ECG

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OV 115 SLIM

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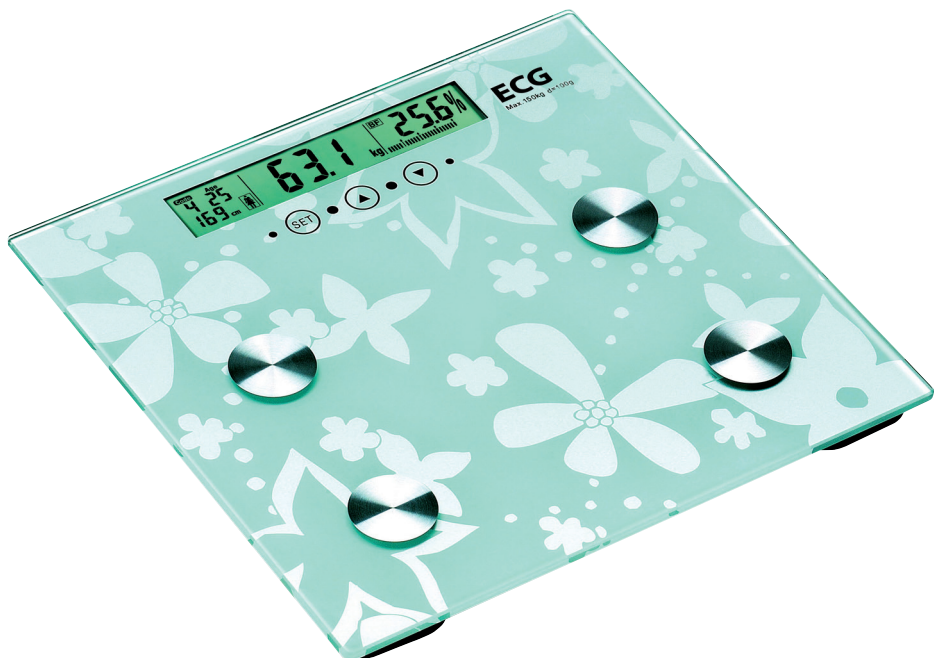
OV 115 SLIM

USER'S MANUAL

Read this user's manual and the safety instructions contained in it carefully before putting the product into operation. The user's manual must always come with the product.

SAFETY INSTRUCTIONS

Read carefully and store for future use!

Warning: The safety instructions in this user's manual do not include all the possible conditions and situations which may arise. The user must understand that a factor that cannot be built into any product is common sense, caution and care. These factors must therefore come from the user/s using and operating this product.

We do not take responsibility for damages caused during transport, by incorrect use or a change or modification of any part of the product.

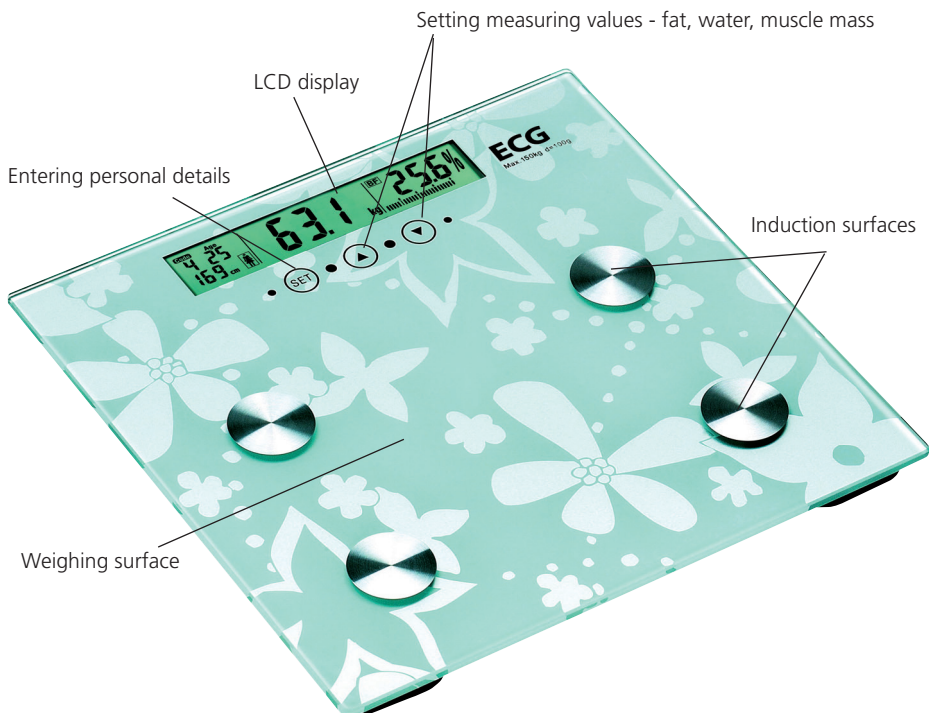
To prevent injury, basic provisions should be adhered to during use, including the following:

1. The personal scale OV 115 SLIM is designed exclusively for household use.
2. **Have all repairs done by a qualified technician!**
3. Do not use the scale in the vicinity of devices with an electromagnetic field. Place the scale sufficiently far away from such a device to prevent the devices from interfering with each other.
4. The result is only an approximate value. This product is not a medical device. Please, always consult with your doctor the medication or diet necessary to achieve optimal weight.
5. Do not place any objects into the openings of the product.
6. The scale is designed for weighing people from 10 years of age.
7. Pets must not be weighed on the scale.
8. Protect the scale from being dropped and impacts and do not jump on it.
9. You or the product may slip on wet surfaces. Likewise on a smooth surface never step on the edge of the scale and while weighing maintain balance on the scale.
10. For reason of electronic impedance **do not use the product in the following cases:**
 - a) if you have an implanted cardio-stimulator (pace maker) or another electronic medical implant
 - b) you have heart problems or disease
 - c) during pregnancy
 - d) if you undergoing dialysis treatment
 - e) you have a fever
 - f) you suffer from osteoporosis
 - g) you have an oedema or swelling
 - h) if you have taken part in exercise lasting longer than 10 hours
 - i) your heart pulse is slower than 60 beats per minute
11. Use the scale only in accordance to the instructions contained in this user's manual.
12. The manufacturer takes no responsibility for damages caused by the incorrect use of this device.
13. This product may not be used by any persons (including children) with mental disabilities, weak strength or by inexperienced persons, unless they have been properly trained or schooled in the safe use of the product, or unless they are properly supervised by a qualified person who is responsible for their safety. Children must be supervised to ensure that they do not play with the product.

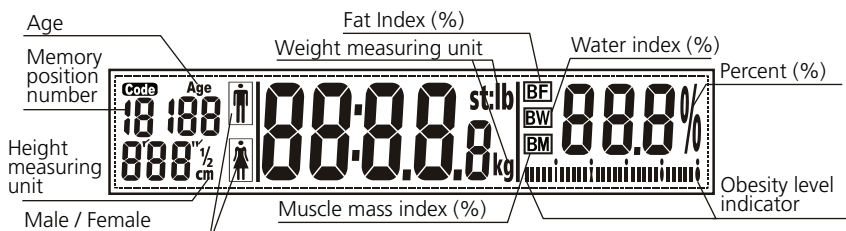
Attention

1. Before weighing always remove your shoes and socks.
2. It is recommended to carry out weighing every day at the same time.
3. Measured results may be misleading after intensive exercise, a significant diet or during extremely high temperatures.
4. Measured values for fat, water and muscle mass may have a significant deviation for the following people:
 - children younger than 10 years of age and adults older than 100 years of age (they can use the device only for basic weighing of body weight),
 - adults over 70 years of age,
 - people operating a fitness business or professional athletes.

DESCRIPTION

LCD display

Touch buttons

-  Turn on the device and select other details (see chapter Setting personal details).
-  Every subsequent time that this button is touched the value of the detail increases by one measuring unit.
-  Every subsequent time that this button is touched the value of the detail decreases by one measuring unit. Note: Holding down the button will speed up the setting of the value of the details.

Changing batteries

1. First remove the insulation foil from the contact of the new battery.
2. Open the battery compartment cover on the underside of the scale.
3. Using a sharp object remove the used battery (e.g. using tweezers).
4. Insert the new battery so that one side of the battery is inserted under the latch and then press on the other side. Make sure to maintain the correct polarity.

After the batteries are inserted for the first time and it is turned on the following details will be displayed:

Gender	Height	Age
Male	165 cm/5'05.0"	25

BASIC MEASUREMENT OF WEIGHT

The scale must be placed on an even, hard and stable surface (not on carpet, the weighing would be inaccurate).

1. Carefully step on the scale - the scale will turn on automatically. Stand evenly on the scale, do not move and wait until your weight is shown on the display has stabilised and locked itself.
2. When you step off the scale it will turn itself off automatically. The device will be automatically turned off when the display shows "0.0" or after about 8 seconds if it is displaying the same weight.

Adding weight of by additional weighing

If you want to find the weight of your child for example, first weigh yourself (wait until the digits have stabilised), then without stepping off the scale pick up your child. The digits will start to flash and then calm down and the shown weight will stabilise. Now you will see the total weight. By subtracting your own weight from the newly identified total value you will get the difference, which represents the weight of your child.

Note: The object being weighed must have more than 2 kg.

MEASURING BODY FAT, WATER AND MUSCLE MASS

Measuring principles

Electric signals go through the human organism at certain frequencies, and the impedance of body fat is higher than the impedance of muscle mass or other human tissues. The device generates an alternating electrical current at a safe frequency and low voltage. Using an analysis of biological resistance, which acts on the passage of this current through the human body (the impedance of various tissues, fat and muscles varies), on the basis of the resistance differences the percentage content of body fat in the human organism is deduced.

To % of body fat content

% of fat is an indicator of the percentage content of body fat in the human organism.

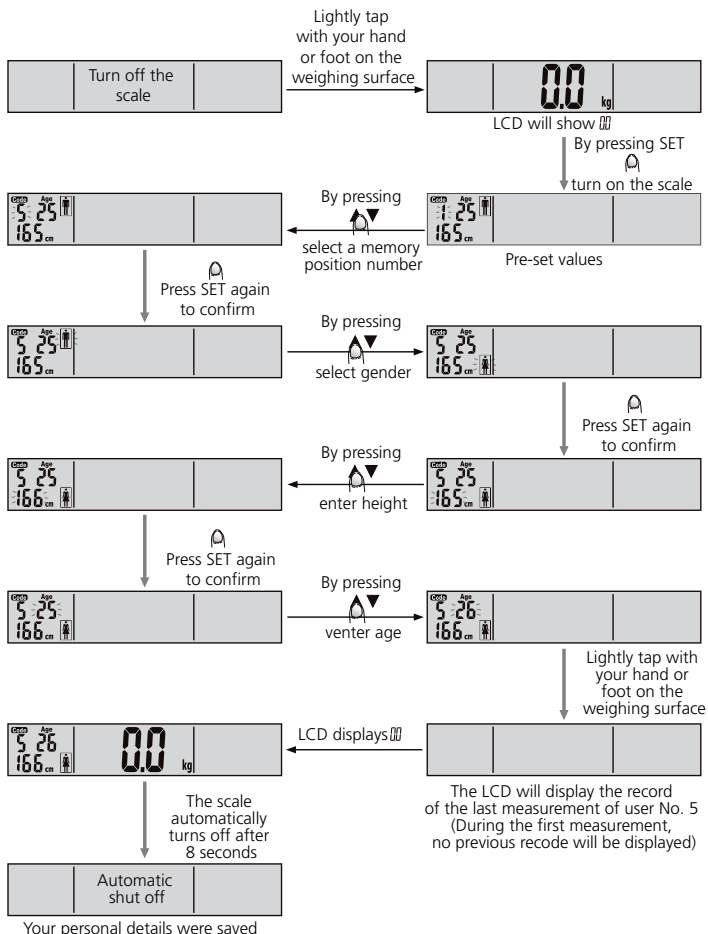
Body fat in the human body should be at a certain level, excessive body fat worsens metabolic processes and increases the risk of various illnesses.

1. Setting personal details

Personal details can be set in the following ranges:

Memory	Gender	Height	Age
1 ~ 12	male/female	100 ~ 250 cm / 3'03.5" ~ 8'02.5"	10 ~ 100

Example: User No. 5, female 166 cm, 26 years old, can set her details as follows:

Note:

- While entering personal details you can gently tap the weighing surface with your hand or foot to set the scale to zero. You can begin measuring body fat, water and muscle mass with the earlier entered personal details or you can update the details originally entered.
- After completing the entry of personal details you can press SET again and continue at any time with entering personal details for another memory position.

2. Measuring body fat, water and muscle mass

Example: User No. 5, female 166 cm, 26 years old, can measure herself as follows:

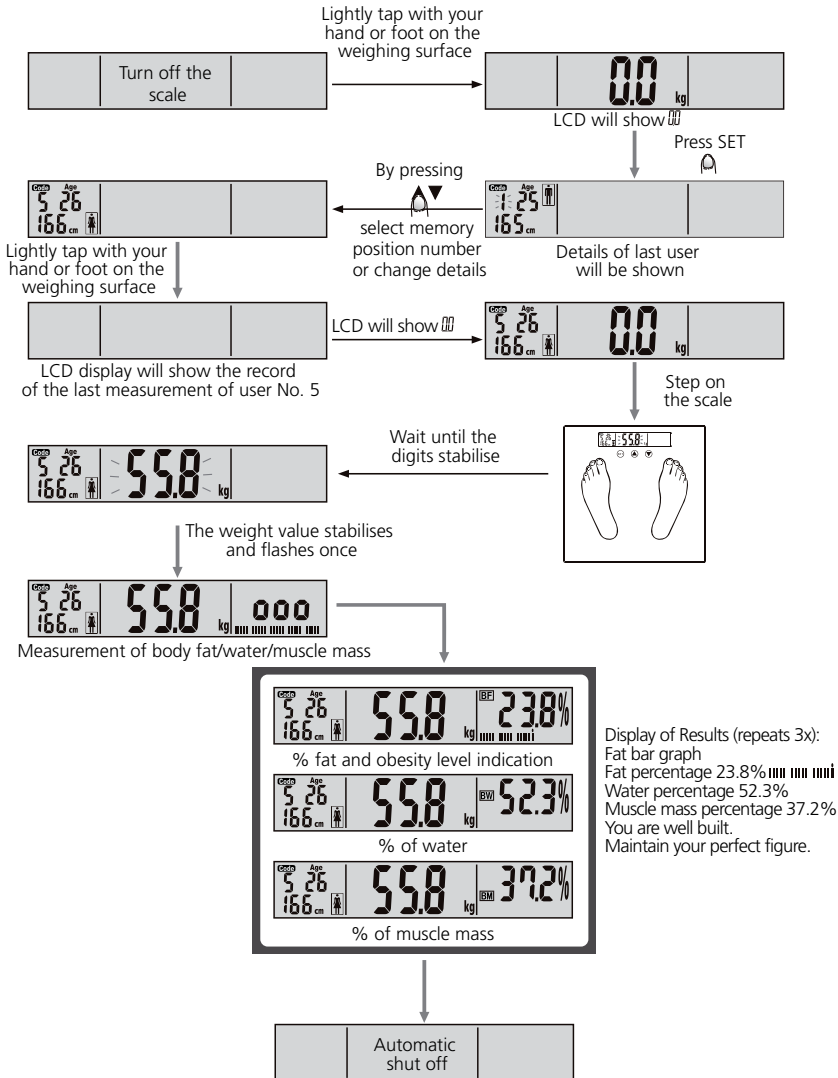
Table of body condition according to fat/water content

Age	Fat (F)	Share of Water(F)	Share of Fat (M)	Share of Water (M)	Condition state	Indication level of obesity
≤30	4.0–16.0 %	66.0–57.8 %	4.0–11.0 %	66.0–61.2 %	too slim	
	16.1–20.5 %	57.7–54.7 %	11.1–15.5 %	61.1–58.1 %	slim	
	20.6–25.0 %	54.6–51.6 %	15.6–20.0 %	58.0–55.0 %	normal	
	25.1–30.5 %	51.5–47.8 %	20.1–24.5 %	54.9–51.9 %	obese	
	30.6–60.0 %	47.7–27.5 %	24.6–60.0 %	51.8–27.5 %	very obese	
>30	4.0–20.0 %	66.0–55.0 %	4.0–15.0 %	66.0–58.4 %	too slim	
	20.1–25.0 %	54.9–51.6 %	15.1–19.5 %	58.3–55.3 %	slim	
	25.1–30.0 %	51.5–48.1 %	19.6–24.0 %	55.2–52.3 %	normal	
	30.1–35.0 %	48.0–44.7 %	24.1–28.5 %	52.2–49.2 %	obese	
	35.1–60.0 %	44.6–27.5 %	28.6–60.0 %	49.1–27.5 %	very obese	

Note: Values in this table are only of informative character.

WARNING MESSAGES

Flat battery

The battery is flat, replace with a new battery.

Overload

Scale overloaded - the maximum load of 150 kg is displayed on the scale.

Please step off the scale to prevent damaging it.

Body fat % is too low

Body fat % is too low. Increase the nutrition in your diet and take good care of yourself.

Body fat % is too high

Body fat % is too high. Watch your diet and exercise more.

Repeated measurement

If after stepping of the scale "C" is displayed, it means that there was some kind of error in the weight measurement.

Please repeat the weighing/measurement.

MAINTENANCE AND CLEANING

1. The glass surface becomes slippery if damp / wet. Keep it dry!
2. Do not move while weighing.
3. Do not subject the scale to impacts, shaking and do not drop it on the floor.
4. The scale is a precise device requiring careful handling.
5. To clean the scale use a damp cloth and make sure that no water gets inside the scale.
DO NOT USE any chemical cleaning agents / abrasive cleaning agents. Do not submerge the scale in water!
6. Use the scale in a cool room with dry air.
7. Always keep the scale in a horizontal position.
8. If you are not successful in turning on the scale, check whether the battery is inserted or whether it is flat. Insert the battery again or replace with a new one.
9. If an error is displayed on the display or if it is not possible to turn off the scale for a long time, please, take out the battery for approximately 3 seconds and then put it back. If you are not successful in solving this problem, please, the store where you bought the device for advice.

TECHNICAL INFORMATION

Fitness scale for measuring body fat, water and muscle mass

Made from 6mm safety glass

LCD display

Touch buttons

Measurement: body weight 2–150 kg
body fat content 4,0–60,0 %
body water content 27,5–66,0 %
Muscle mass content 20–56 %

Weighing sensitivity 100 g

Fat and water measurement sensitivity 0.1%

Memory for 12 people

Automatic zeroing

Automatic on and off

Indication of overload and weak batteries

Power source: 2 x 3V lithium batteries (included)

USE AND DISPOSAL OF WASTE

Wrapping paper and cardboard - take to waste collection point. Wrapping foil, PE bags, plastic parts - take to plastic recycling bins. Hand over flat batteries to the retailer or take to a place defined by the municipality for the disposal of dangerous wastes.

08/05

DISPOSAL OF PRODUCTS AT END OF LIFETIME

The disposal of used electric and electronic devices (valid in EU member states and other European countries with an implemented system for sorting waste).

The displayed symbol on the product or the packaging means that the products should not be handled as household waste. Hand over the product to a location designated for the recycling of electric and electronic devices. By correctly disposing of the product you will prevent negative effects on human health and the environment. The recycling of material helps to protect natural resources. You can receive more information about the recycling of this product at the town council, at the organisation for the processing of household waste or at the store where you purchased this product.

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© ZStudio 2010
K+B E-Tech GmbH & Co. KG